



# ColorFlex<sup>®</sup> EZ



measure  
color  
the way your  
eye sees it



## DESIGNED FOR

Coffee  
powders  
GRANULES  
Sugar  
PELLETS  
Liquids  
detergents  
resins  
pigments  
CHIPS  
laminates  
sauces  
grout  
Soups  
SPICES  
film  
cement &  
concrete  
YOGURTS  
juices  
panels &  
plaques  
adhesives  
COLORANTS  
all dry chemicals  
DAIRY  
paint panels  
fruit  
Cheeses  
powder  
coatings  
nuts  
SOAPS  
Solids

## The Power to See Color the Way Your Customers Do

HunterLab's next generation ColorFlex® EZ spectrophotometer takes color quality control to its highest level with 45°/0° design for the ultimate in color measurement precision.

Combining versatility, simplicity and performance, the ColorFlex® EZ reflects

### over 60 years of color measurement innovation

in one easy-to-use, compact instrument from the world's most trusted experts in color quality.

## 45°/0° Design: For Relentless Perfection in Color Quality

HunterLab's ColorFlex® EZ offers the advantage of 45°/0° design ensuring that you not only get your colors right - but perfect every time. It is the optimal and only instrument design that measures your samples the way the human eye does. It lets you see your colors exactly the way your customers do, not just in the lab but in the real world.

## More Connectivity


More options with 3 USB 2.0 communications ports to connect to a wide variety of devices:

- USB flash drives to transfer results to several machines or computers
- Keyboards to enter custom sample IDs manually
- Barcode scanners to acquire samples fast without data entry errors
- Printers to get instant hard copies of sample readouts
- Computers to plot data using HunterLab's optional EasyMatch® QC software

## More Simplicity

Easy to operate for all your color measurement challenges:

- Small size requiring very little benchtop space
- "One-touch" operation

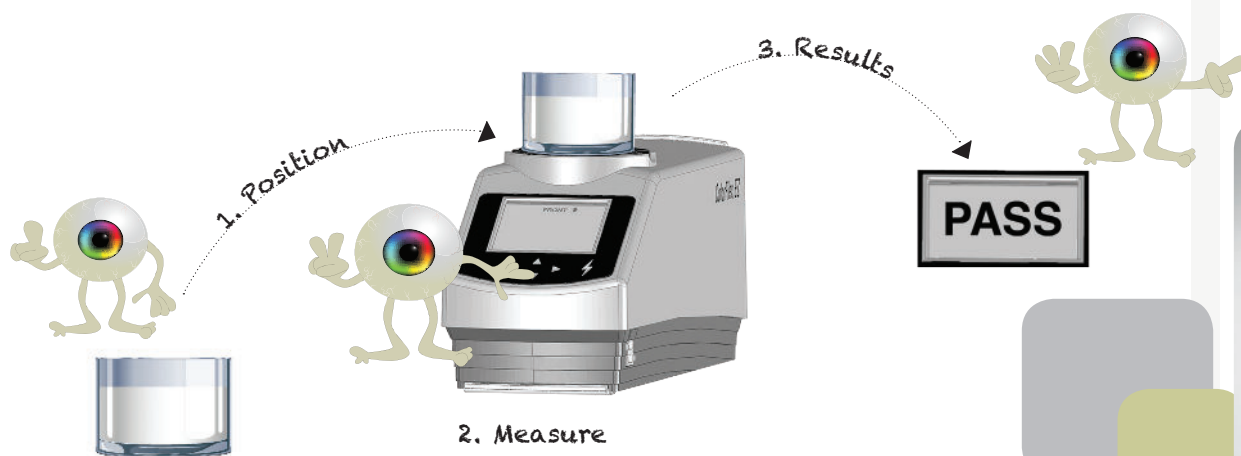

## More Performance

More capabilities for increased performance and better outcomes:

- Sampling versatility to measure everything from opaque solids, liquids, powders, granules and pellets to translucent solids and liquids
- Easy data transfer to move set ups and standards from one instrument to another and upload or download to USB flash drives for back ups
- **Expanded sample storage for 250 setups, product standards and tolerances**  
and up to **2000 sample measurements**



## More Flexibility

Built-in software using 250 set up configurations with product standard values and tolerances for a wide variety of parameters including:

- Color Scale, Indices, Illuminant and Observer to specify conditions under which data is reported
- Standard values to identify product standards or targets for color comparisons
- Pass / fail tolerances to measure against product specifications
- Averaging to provide accurate measurements of uneven or inconsistent samples
- Display formats to report measurement data multiple ways such as color data, color difference data, and spectral difference data
- A variety of sample handling accessories

Also available  
ColorFlex® EZ Tomato

## The HunterLab Advantage

The ColorFlex® EZ is backed by over 60 years of quality innovation and experience from HunterLab, the world's most trusted color quality experts.


With an unmatched reputation for delivering the right solution for any challenge, HunterLab tailors products and technologies for every color measurement need and budget offering the broadest range of color measurement solutions in the industry.

It's that easy!  
It's ColorFlex® EZ



For more information go to [www.hunterlab.com](http://www.hunterlab.com)

or contact your local HunterLab representative.

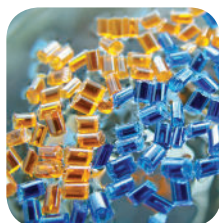

# Accessories



04-4000-00

## Sample Cup Opaque Cover

Provides a light trap to exclude the interference of external light on the sample cup.

04-7209-00

## Glass Sample Cup: 64 mm (2.5 in)

Optically clear glass cup for sample presentation of liquids, powders, granules and pellets.

02-4579-00

## Ring and Disc Set

For use in the glass sample cup for sample presentation of liquids, powders, granules and pellets.



04-6622-00

## Port Insert for Sample Cup

**64 mm (2.5 in)** Designed to hold and position the glass sample cup perfectly to accurately measure powders, granules, pellets, paste and liquids.

04-6623-10

**32 mm (1.25 in)** Port Insert

04-6623-09

**25 mm (1 in)** Port Insert

04-6623-08

**19 mm (0.75 in)** Port Insert

04-6623-05

**13 mm (0.50 in)** Port Insert

Enable you to adjust the sample measurement area.



A13-1014-294

## USB Flexible Keyboard

88 character keyboard allows user to enter IDs directly into ColorFlex® EZ.

A13-1014-254

## USB Barcode Reader

Barcode Scanner scans product IDs directly into ColorFlex® EZ. Automatically detected at USB port.

A13-1014-259

## USB Printer

USB printer allows hard copy measurements to be printed. Requires a standard USB cable, sold separately. Uses Thermal Printer paper, sold separately.



L02-1014-608

## Hard-Sided Carrying Case, ColorFlex® EZ

Hard-sided Carrying Case for ColorFlex® EZ.